

Rapid Response™

Clostridium difficile GDH Test Cassette (Feces)

REF CdifAB-9E25

Product Insert

In vitro rapid diagnostic test for the detection of Clostridium difficile GDH antigen in human feces samples.

For professional *in vitro* diagnostic use only.

Intended Use

The Rapid Response™ *Clostridium difficile* GDH Test Cassette is a rapid chromatographic immunoassay for the qualitative detection of *Clostridium difficile* GDH antigen in the human feces specimen.

Summary

Clostridium difficile is an anaerobic bacterium acting as an opportunistic pathogen: it grows in the intestine when the normal flora has been altered by treatment with antibiotics.^{1,2,3} Toxinogenic strains of *Clostridium difficile* cause infections from mild-diarrhea to pseudomembranous colitis, potentially leading to death.⁴

Disease is caused by two toxins produced by toxinogenic strains of *C. difficile*: Toxin A (tissue-damaging enterotoxin) and Toxin B (cytotoxin). Some strains produce both toxins A and B, some others produce Toxin B only. The potential role of a third (binary) toxin in pathogenicity is still debated.⁴

The use of Glutamate Dehydrogenase (GDH) as an antigen marker of *C. difficile* proliferation has been shown to be very effective because all strains produce high amount of this enzyme.^{5,6}

The Rapid Response™ *Clostridium difficile* GDH Test Cassette allows the specific detection of *C. difficile's* GDH in stool specimen. Samples with a positive result should be investigated further to test for toxigenicity of the bacteria.

Principle

This is a ready-to-use test that is based on the use of a membrane technology with colloidal gold. A nitrocellulose membrane is sensitized with antibody directed against *Clostridium difficile* antigen (GDH). The test's specificity is ensured by an antibody specific to the *Clostridium difficile* GDH that is conjugated to the colloidal gold. This conjugate is dried on polyester.

The fecal sample must be diluted into the extraction buffer that is supplied with the test. When the fecal suspension come into contact with the strip, the solubilized conjugate migrates with the sample by passive diffusion and the conjugate and sample material come into contact with the anti-Clostridium antibody

adsorbed on to the nitrocellulose. If the sample contains *C. difficile* GDH, the conjugate-antigen complex will remain bound to the anti-*C. difficile* GDH reagent and a red line will develop. Solution continues to migrate to encounter a second reagent that binds the migration control conjugate, thereby producing a red control line that confirms that the test is working properly. The result is visible within 10 minutes.

Reagents

The test cassette contains anti-*Clostridium difficile* GDH particles and anti-*Clostridium difficile* GDH coated on the membrane.

Precautions

- All operations linked to the use of the test must be performed in accordance with Good Laboratory Practices (GLP).
- All reagents are for *in vitro* diagnostic use only.
- Avoid touching nitrocellulose with your fingers.
- Wear gloves when handling samples.
- Never use reagents from another kit.
- Reagents' quality cannot be guaranteed beyond their shelf-life dates or if reagents are not stored under required conditions as indicated in the insert.
- Dispose of gloves, swabs, test tubes and used devices in accordance with GLP.
- Each user is responsible for the management of any waste produced, and must ensure that it is disposed of in accordance with the applicable legislation.

Materials

Materials provided

- Test cassettes
- Specimen collection tubes with buffer
- Droppers
- Product insert

Materials required but not provided

- Timer
- Specimen collection container

Storage and Stability

Store as packaged at room temperature or refrigerated (35.6-86°F; 2-30°C). The test is stable through the expiration date printed on the sealed pouch. The test must remain in the sealed pouch until use. **DO NOT FREEZE.** Do not use beyond the expiration date.

Collection and Storage of Specimens

The stool specimens must be tested as soon as possible after collection. If necessary, original feces specimen could be stored at 35.6-46.4°F (2-8°C) for 3 days or -4°F (-20°C) for longer periods of time; extracted specimen in buffer could be stored at 35.6-46.4°F (2-8°C) for 1 week or -4°F (-20°C) for longer

periods of time.

Make sure that the specimens are not treated with solutions containing formaldehyde or its derivatives.

Test Procedure


Allow the test, specimen, stool collection buffer and/or control to equilibrate to room temperature (59-86°F; 15-30°C) prior to testing.

1. To process fecal specimens:

- For **Solid Specimens:**
 - Unscrew the cap of the specimen collection tube, then **insert the specimen collection applicator into the fecal specimen at least 3 different sites** to collect approximately **50 mg of feces** (equivalent to 1/4 of a pea). Do not scoop the fecal specimen.
- For **Liquid Specimens:**
 - Hold the dropper vertically, aspirate fecal specimens, and then transfer **2 drops of the liquid specimen** (approximately 80 µL) into the specimen collection tube containing the extraction buffer.

Tighten the cap onto the specimen collection tube, then **shake the specimen collection tube vigorously** to mix the specimen and the extraction buffer. Leave the collection tube for reaction for 2 minutes.


- Remove the test cassette from the foil pouch and use it as soon as possible. Best results will be obtained if the test is performed immediately after opening the foil pouch.
- Hold the specimen collection tube upright and **unscrew the cap** of the specimen collection tube. Invert the specimen collection tube and **transfer 3 full drops of the extracted specimen** (approximately 120µL) to the specimen well (S) of the test cassette, then start the timer. Avoid trapping air bubbles in the specimen well (S). See illustration below.
- Read the results at **10 minutes** after dispensing the specimen. Do not read results after 20 minutes.
- NOTE:** If the specimen does not migrate (presence of particles), centrifuge the diluted sample contained in the extraction buffer vial. Collect 80 µL of supernatant, dispense into the specimen well (S) of a new test cassette and start afresh following the instructions mentioned above.

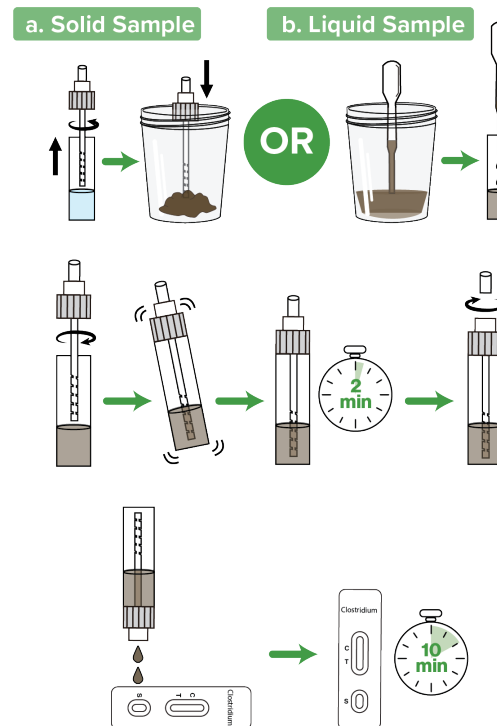


Results Interpretation


POSITIVE: * Two colored lines appear. One colored line should be in the control line region (C) and another apparent colored line should be in the test line region (T).

***NOTE:** The intensity of the color in the test line region (T) will vary depending on the concentration of Clostridium difficile GDH antigen present in the specimen. Therefore, any shade of color in the test line region (T) should be considered positive.

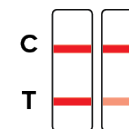

NEGATIVE: One colored line appears in the control line region (C). No line appears in the test line region (T).

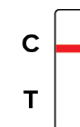

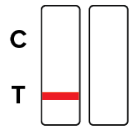

INVALID: Control line (C) fails to appear. Insufficient specimen volume or incorrect procedural techniques are the most likely reasons for control line failure. Review the procedure and repeat the test with a new test cassette. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

NOTE: During the drying process, a very faint shadow may appear at the Test line position. It should not be regarded as a positive result.

Quality Control

An internal procedural control is included in the test. A colored line appearing in the control line region (C) is an internal positive procedural control. It confirms sufficient specimen volume, adequate membrane wicking and correct procedural technique. Control standards are not supplied with this kit; however, it is recommended that positive and negative controls be tested as a good laboratory practice to confirm the test procedure and to verify proper test performance.

Limitations

1. The test is qualitative and cannot predict the quantity of antigens present in the sample. Clinical presentation and other test results must be taken into consideration to establish diagnosis.
2. A positive test does not rule out the possibility that other pathogens may be present.

Performance Characteristics

Detection Limit

The detection limit was evaluated by diluting a purified GDH preparation and the results show that the concentration of protein detected is 1ng/mL.

Sensitivity – Specificity

Method	Results	Other Rapid Test		Total Results
		Positive	Negative	
Rapid Response™	Results			
<i>Clostridium difficile</i>	Positive	78	2	80
GDH Test Cassette	Negative	1	119	120
Total Results		79	121	200

Relative Sensitivity: 98.7% (95%CI:*93.1%-100%)

Relative Specificity: 98.3% (95%CI:*94.2%-99.8%)

Relative Accuracy: 98.5% (95%CI:*95.7%-99.7%)

*Confidence Intervals

Repeatability and Reproducibility

To check intra-batch accuracy (repeatability), the same positive samples and a buffer solution were processed 15 times on kits of the same production batch in the same experimental conditions. All observed results were confirmed as expected. To check inter-batch accuracy (reproducibility), some samples (positive and buffer) were processed on kits from three different production batches. All results were confirmed as expected.

Cross Reactivity


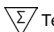






An evaluation was performed to determine the cross reactivity of the Rapid Response™ *Clostridium difficile* GDH Test Cassette. No cross reactivity against gastrointestinal pathogens occasionally present as following:

Campylobacter coli	Salmonella enteritidis
Shigella dysenteriae	Campylobacter jejuni
Salmonella paratyphi	Shigella flexneri
E.coli O157:H7	Salmonella typhi
Shigella sonnei	H.pylori
Salmonella typhimurium	Staphylococcus aureus
Listeria monocytogenes	Shigella boydii
Yersinia enterocolitica	

Bibliography

1. RamadassBalamurugan, V. Balaji and Balakrishnan S. Ramakrishna: Estimation of faecal carriage of *Clostridium difficile* in patients with ulcerative colitis using real time polymerase chain reaction, Indian Journal of Medical Research, p.472-477, May 2008
2. E. J. Kuijper, B. Coignard and P. Tüll: Emergence of *Clostridium difficile*-associated disease in North America and Europe, Review Clinical Microbiology and Infections, 12 suppl6, p. 2-18, Oct. 2006
3. Leyerly D.M., H.C. Krivan and D.T. Wilkins: *Clostridium difficile*: its disease and toxins. Clinical Microbiology Reviews, p. 1-18, Jan. 1988
4. Ramsey L. et al: Fulminant *Clostridium difficile*: an underappreciated and increasing cause of death and complications, Annals of Surgery 235 (3) p. 363-372: Mar. 2002
5. Wren MW., Kinson R., Sivapalan M., Shemko M., Shetty NR.: Detection of *Clostridium difficile* infection: a suggested laboratory diagnostic algorithm, British Journal of Biomedical Sciences, 66(4) p. 175-179, 2009.
6. Willis DH. And JA Kraft: Confirmation that the latex-reactive protein of *Clostridium difficile* is a Glutamate Dehydrogenase. Journal of clinical microbiology, 30, p. 1363-1364, May 1992

Glossary of Symbols

	Consult instructions for use		Test per Kit		Do Not Reuse
	Store between 35.6°F to 86°F		Use by		Catalogue #
	Lot Number		For <i>in vitro</i> diagnostic use only		

BTNX Inc.
722 Rosebank Road,
Pickering, ON L1W 4B2
Canada



Technical Support: 1-888-339-9964