

CLIA Waived Rapid Response™ Urine Drug Cups

DRUG SCREENING

CPT Code: 80305QW

An FDA-cleared, CLIA-Waived one-step solution for urine drug screening.

Rapid Response™ Drug Test Cups simultaneously detect multiple drugs in urine samples. With a rapid 5-minute turnaround time, our multi-drug cups empower your team to quickly identify patient needs and initiate counseling during the same visit, enhancing patient engagement. Available with a wide range of parameters and two user-friendly formats: the traditional Round Cup and the Flat Cup.

Sample: Urine

Time-to-Result: 5 minutes

Storage Condition: 35.6-86°F / 2-30°C

Test Principle: Lateral Flow Assay

For *in vitro* diagnostic use only

Visit our Website

Custom pricing and combinations available for bulk orders, contact info@lochnessmedical.com for details.



Detects Multiple Drugs

This CLIA-Waived drug screening cup detects a wide range of drugs in various combinations, including fentanyl. It's an ideal choice for professionals seeking a precise and reliable drug screening solution.



Point-of-Care Testing

Screen for multiple drugs and metabolites in patient samples in just 5 minutes.



Easy & Convenient

A simple one-step test with up to 16 integrated drug test strips enables easy collection and testing. Each cup features a leak-proof cap and a temperature strip for added reliability.



Kit Content


Each kit includes test cups, an adulteration color chart (when applicable), and a product insert.

Test Procedure

1. Open the cap and add specimen. Specimen volume must exceed the minimum line.
2. While the cup is on a flat surface, replace and secure the cap on the cup.
3. 4 minutes after specimen collection, check the temperature label.
4. Peel off the label on the test cup to read results at 5 minutes.

Test procedure is the same for both the Round Cup and Flat Cup formats.

These instructions are for illustrative purposes only. Please read the instructions supplied with the test before use.

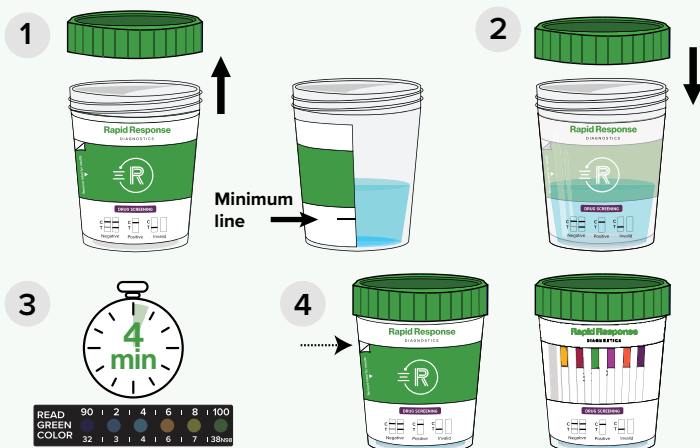

**Custom Combinations available,
Contact us for details.**



Round Cup



Traditional round urine drug test cup featuring a leak-proof cap, a simple procedure, and results in just 5 minutes. Fentanyl testing is now available with select combinations.

Round Cups	Product Code	Drug Parameters
Multi-Drug One Step Cup - 10 Drugs	D10.14-1V	COC300 AMP1000 MET1000 THC50 MTD300 OPI2000 OXY100 PCP25 BAR300 BZO300
Multi-Drug One Step Cup - 12 Drugs with Fentanyl and Adulterant Tests	D12.49-1VA	FYL1 OXY100 BZO300 OPI300 MET500 AMP500 COC150 THC50 MDMA500 BUP10 BAR300 MTD300 + Oxi/Cre/pH
Multi-Drug One Step Cup - 13 Drugs with Fentanyl Test	D13.5-1V	FYL1 OPI300 MDMA500 OXY100 PCP25 MTD300 BAR300 BZO300 COC150 AMP500 MET500 THC50 BUP10

Flat Cup



Same simple procedure and fast turnaround as the Round Cup, but compatible with the RapidReader App and Lateral Flow Scanner. Read more below.

Flat Cups	Product Code	Drug Parameters
Multi-Drug One Step (Flat) Cup - 5 Drugs	D5.1-1VF	COC300 AMP1000 MET1000 THC50 OPI2000
Multi-Drug One Step (Flat) Cup - 5 Drugs	D5.75-1VF	COC300, MET1000, AMP1000, OPI2000, OXY100
Multi-Drug One Step (Flat) Cup - 12 Drugs	D10.3-1VF	COC300, AMP1000, MET1000, THC50, OPI2000, PCP25, BAR300, BZO300, MTD300, MDMA500
Multi-Drug One Step (Flat) Cup - 10 Drugs	D10.32-1VF	COC300, MET1000, AMP1000, OPI2000, OXY100, MTD300, PCP25, BAR300, BZO300, MDMA500
Multi-Drug One Step (Flat) Cup - 12 Drugs with Adulterant Tests	D12.11-1VAF	MOP300, OXY100, BZO300, OPI2000, MET1000, AMP1000, COC300, THC50, MDMA500, BUP10, BAR300, MTD300, + Cre/pH/Oxi
Multi-Drug One Step (Flat) Cup - 12 Drugs	D12.5-1VF	AMP1000 BAR300 BUP10 BZO300 COC300 MDMA500 MTD300 MET1000 OPI2000 OXY100 THC50 PCP25
Multi-Drug One Step (Flat) Cup - 14 Drugs	D14.5-1VF	AMP1000 BAR300 BUP10 BZO300 COC300 MET1000 MDMA500 OPI2000 MTD300 OXY100 PCP25 THC50 PPX300 TCA1000

 Compatible with
The RapidReader App

 Compatible with the
Lateral Flow Scanner

The RapidReader App



The RapidReader App helps treatment centers streamline their drug screening process through a compact and user-friendly mobile application. Using a cellphone or tablet, the app utilizes AI to capture individual test results from a photo and automatically uploads the report to the patient's profile in the EMR.

Scan or [click here](#)
to learn more.



Lateral Flow Scanner



The Rapid Response™ Lateral Flow Scanner is an ideal solution for treatment centers looking to streamline medium- to high-volume drug screen reporting. The system can simultaneously scan up to three patients' Rapid Response™ Drug Tests. Powered by AI, it captures the results and uploads a report directly to each patient's profile in the EMR system.

Scan or [click here](#)
to learn more.

