


# Rapid Response™ hCG Pregnancy Tests

## WOMENS HEALTH

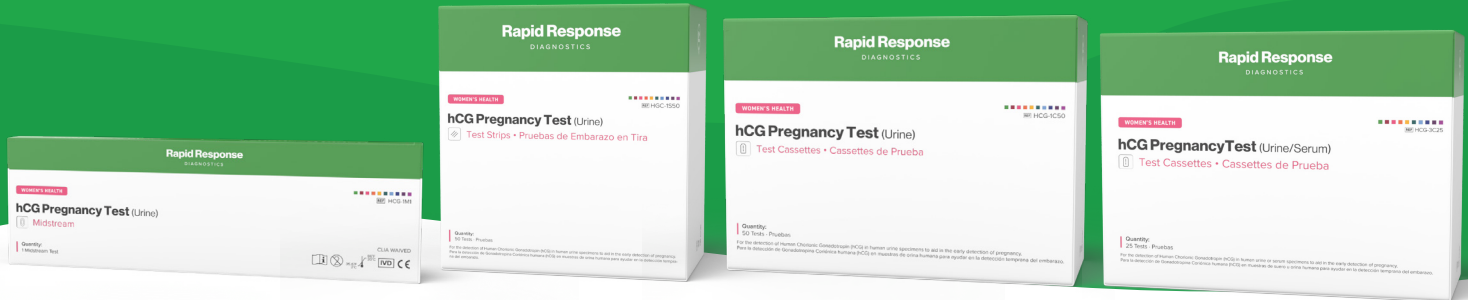

**HCG-1M1**  
hCG Pregnancy Test  
Midstream (Urine)

**Sample:** Urine  
**Format:** Midstream  
**Time-to-Result:**  
 3 minutes  
**Storage Condition:**  
 36-86°F/2-30°C  
**Test Principle:**  
 Lateral Flow Assay  
**CLIA Status:**  
 Waived



**HCG-1S50,  
HCG-1S100**  
hCG Pregnancy Test  
Strip (Urine)

**Sample:** Urine  
**Format:** Strip  
**Time-to-Result:**  
 5 minutes  
**Storage Condition:**  
 36-86°F/2-30°C  
**Test Principle:**  
 Lateral Flow Assay  
**CLIA Status:**  
 Waived

*Available in 50 or 100 pack.*



**HCG-1C25,  
HCG-1C50**  
hCG Pregnancy Test  
Cassette (Urine)

**Sample:** Urine  
**Format:** Cassette  
**Time-to-Result:**  
 3 minutes  
**Storage Condition:**  
 36-86°F/2-30°C  
**Test Principle:**  
 Lateral Flow Assay  
**CLIA Status:**  
 Waived

*Available in 25 or 50 pack.*



**HCG-3C25,  
HCG-3C50**  
hCG Pregnancy Test  
Cassette (Urine/Serum)

**Sample:** Urine/Serum  
**Format:** Cassette  
**Time-to-Result:**  
 5 minutes  
**Storage Condition:**  
 36-86°F/2-30°C  
**Test Principle:**  
 Lateral Flow Assay  
**CLIA Status:**  
 Waived

*Available in 25 or 50 pack.*

With a range of accurate, easy-to-use, early-detection pregnancy tests, the Rapid Response™ Women's Health product line offers medical professional with solutions suited to the unique healthcare and diagnostic needs of women.

The Rapid Response™ hCG Pregnancy Tests provides visual immunoassays for qualitative detection of human chorionic gonadotropin (hCG) in human urine or serum specimens. The Rapid Response™ hCG Pregnancy Tests are capable of detecting pregnancy as early as 1 day after the first missed menses.

**REF HCG-1M1**

# hCG Midstream Test

## Kit Content

- 1 Test
- Individually packed test midstream
- 1 package insert



## Product Details

**Sample:** Urine

**Format:** Midstream

**Time-to-Result:** 3 minutes

**Storage Condition:**


36-86°F/2-30°C

**Test Principle:** Lateral Flow Assay

## Features

- User-friendly design
- Does not require a urine cup
- Built-in control line
- Double antibody-sandwich technique
- Detects at a concentration of 20 mIU/mL or greater


## Testing Procedure

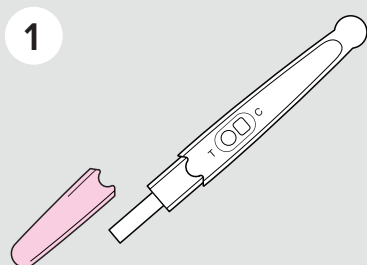

Remove protection cap.



Add sample to the device.

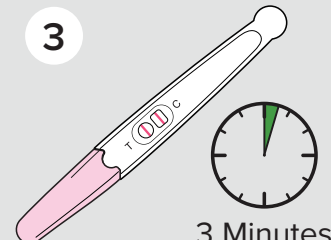

Re-cap protection cap and read results at 3 minutes.

These instructions are for illustrative purposes only. Please read the instructions supplied with the test before use.

## Result Interpretation

### Positive/Pregnant



Two colored bands appear on the membrane. One band appears in the control (C) and another appears in the test region (T).

### Negative/Not Pregnant



Only one colored band appears in the control region (C). No apparent colored band appears in the test region (T).

### Invalid



No line appears in the control (C) region. The result is invalid.

## Ordering Information

Product Code	Product Name	Contents
HCG-1M1	hCG Pregnancy Test	1 Tests


To learn more contact your local **Sales Representative**, call us at **+1 888 506-2658**, or email us at **info@lochnessmedical.com**

REF HCG-1S50, HCG-1S100

# hCG Test Strip

## Kit Content

- 50/100 Test
- Individually packed test strips
- Package insert

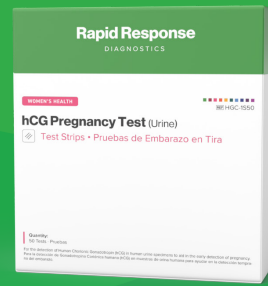

## Product Details

**Sample:** Urine

**Format:** Strip

**Time-to-Result:** 5 minutes

**Storage Condition:**

36-86°F/2-30°C

**Test Principle:** Lateral Flow Assay

## Features

- Compact format makes it easy to transport large numbers for tests
- Internal control line included
- Detects at a concentration of 20 mIU/mL or greater
- Sample can be stored at 36-86°F/2-30°C


## Testing Procedure

1



Remove test from sealed pouch.


2

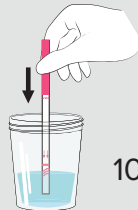

10-15 Seconds

Dip test strip in urine specimen for 10-15 seconds.

3

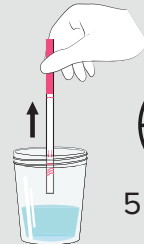

5 Minutes

Remove strip from specimen. Read results at 5 minutes.

These instructions are for illustrative purposes only. Please read the instructions supplied with the test before use.

## Result Interpretation


**Positive**



**Negative**



**Invalid**

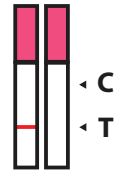

**Accuracy**

**99.8%**

**Relative Sensitivity**

(98.7%-99.9%)\*

**99.9%**

**Relative Specificity**

(99.8%-100.0%)\*

\*95% Confidence Interval

## Ordering Information

Product Code	Product Name	Contents
HCG-1S50, HCG-1S100	hCG Pregnancy Test	50/100 Tests

To learn more contact your local **Sales Representative**, call us at **+1 888 506-2658**, or email us at **info@lochnessmedical.com**

REF HCG-1C25, HCG-1C50

# hCG Test Cassette

## Kit Content

- 25/50 Test
- Individually packed test devices
- Package insert



## Product Details

**Sample:** Urine

**Format:** Cassette

**Time-to-Result:** 3 minutes

**Storage Condition:**

36-86°F/2-30°C


**Test Principle:** Lateral Flow Assay

## Features

- Compact format makes it easy to transport large numbers for tests
- Dual specimen format available
- Internal control line included
- Detects at a concentration of 20 mIU/mL or greater
- Sample can be stored at 36-86°F/2-30°C


## Testing Procedure

1

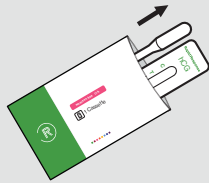


Remove test from sealed pouch.

2

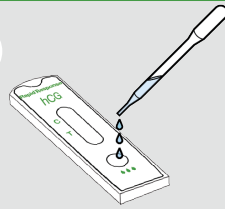

Add 3 drops of specimen into well of cassette.

3

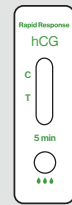


3 Minutes

Read results at 3 minutes.


These instructions are for illustrative purposes only. Please read the instructions supplied with the test before use.

## Result Interpretation


**Positive**

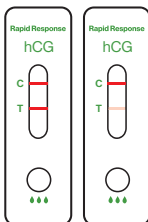

**Negative**

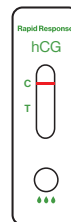

**Invalid**

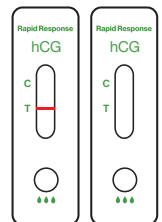

**Accuracy**

**99.8%**

**Relative Sensitivity**  
(98.7%-99.9%)\*

**99.9%**

**Relative Specificity**  
(99.8%-100.0%)\*

\*95% Confidence Interval

## Ordering Information

Product Code	Product Name	Contents
HCG-1C25, HCG-1C50	hCG Pregnancy Test	25/50 Tests

To learn more contact your local **Sales Representative**, call us at **+1 888 506-2658**, or email us at **info@lochnessmedical.com**

REF HCG-3C25, HCG-3C50

# hCG Test Cassette

## Kit Content

- 25/50 Test
- Individually packed test devices
- Package insert



## Product Details

**Sample:** Urine/Serum

**Format:** Cassette

**Time-to-Result:** 5 minutes

**Storage Condition:**

36-86°F/2-30°C


**Test Principle:** Lateral Flow Assay

## Features

- Compact format makes it easy to transport large numbers for tests
- Dual specimen format available
- Internal control line included
- Detects at a concentration of 20 mIU/mL or greater
- Sample can be stored at 36-86°F/2-30°C


## Testing Procedure

1

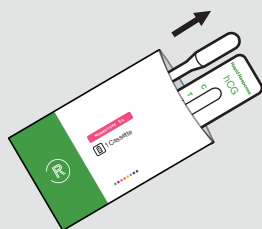


Remove test from sealed pouch.

2

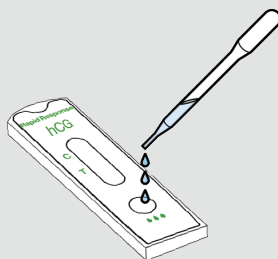

Add 3 drops of urine or serum specimen into well of cassette.

3

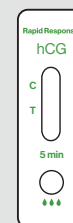


5 Minutes

Read results at 5 minutes.


These instructions are for illustrative purposes only. Please read the instructions supplied with the test before use.

## Result Interpretation


### Positive/Non-Negative

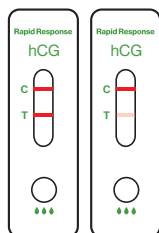

### Negative

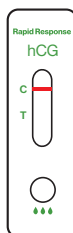

### Invalid

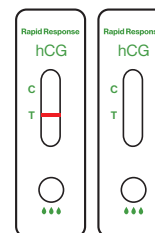

## Ordering Information

Product Code	Product Name	Contents
HCG-3C25, HCG-3C50	hCG Pregnancy Test	25/50 Tests


To learn more contact your local **Sales Representative**, call us at **+1 888 506-2658**, or email us at **info@lochnessmedical.com**

REF HCG-3C25, HCG-3C50


# hCG Test Cassette

## Accuracy

### Urine

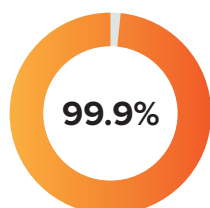


**Relative Sensitivity**  
(97.2%-100.0%)\*

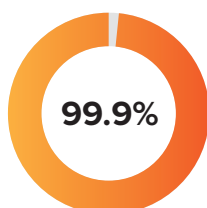


**Relative Specificity**  
(98.0%-100.0%)\*

\*95% Confidence Interval

### Serum

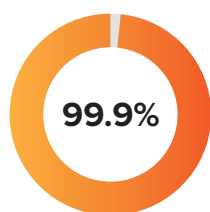

**Relative Sensitivity**  
(97.8%-100.0%)\*

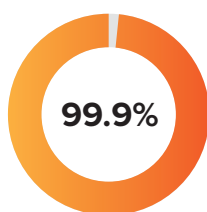

**Relative Specificity**  
(98.5%-100.0%)\*

\*95% Confidence Interval



## Ordering Information

Product Code	Product Name	Contents
HCG-3C25, HCG-3C50	hCG Pregnancy Test	25/50 Tests

To learn more contact your local **Sales Representative**, call us at **+1 888 506-2658**, or email us at **info@lochnessmedical.com**